



## **RELIABLE - QUIET - EASY TO USE**

# Quattro

Bedpan Disposal System

With proven reliability, easy to clean surfaces, classleading quietness and efficiency, the **Quattro** by Haigh has been designed with Nurses, Estate Managers and Site Operators front of mind



## Reliable. For lower running costs.

We absolutely understand that your bedpan machine must work without fail tens, if not hundreds of times per day in a busy sluice-room.

**Quattro** is tested on pulp disposable items from all of the many manufacturers and types of product, from urine bottles to large wash-bowls and sample pots

Haigh test all machines rigorously as they leave our production line. Many features have been expertly built in specifically to reduce the risk of machine failure or blockage of your drainage system:

- ✓ Weir technology controls particle size
- ✓ Automatic and mechanical operation
- ✓ Ease of access on all key components

Foreign objects (gloves etc.) are held back in the **Quattro**, macerated pulp particle size is also regulated by our unique concentric weir technology, controlling particle size and protecting drains from blockage.



**SUPERFINE MACERATION** 

**PROTECT YOUR DRAINS** 

#### QUALITY ASSURED ON ALL TYPES OF PULP

#### Compliant. Best practice.

Your sluice room may not be the most glamorous part of your site but it need not be unpleasant to work in and needs to enable compliant with relevant regulations.

Quattro improves your site environment:

- ✓ Easy to clean
- ✓ Embedded antimicrobial silver ions
- ✓ Efficient use of Energy and Water for ISO14001
- ✓ Quality lid seal restricts odours
- ✓ End of cycle sanitization
- ✓ WRAS approved

Cleaner surfaces combined with the single use pulp concept demonstrably leads to lower risk of cross infection.






## Easy to Use

**Quattro** is designed to make the life of the user as straightforward and pleasant as possible. As the 5<sup>th</sup> generation of the Haigh Sluicemaster range **Quattro** does this by giving you:


- ✓ The lowest loading height
- ✓ The deepest opening
- ✓ High visibility foot pedal opening
- ✓ Automatic starting
- ✓ Fast cycle time
- ✓ The quietest macerator as tested by Quietmark

As a result, not only does **Quattro** aim make the life of the carer better but particularly with its class-leading low sound levels and short run time, the disposal of patient waste is kept discrete with minimal disruption, even during the night.

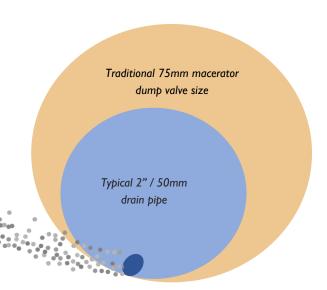


## **Quattro: THE EFFICIENT CHOICE**

### Minimise Water Use

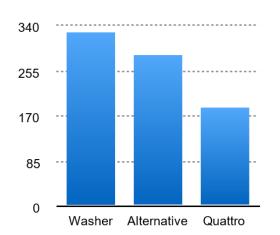

## Minimise Energy Use

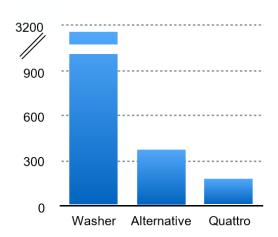

Quattro allows water use to be optimised for your installation.

- Minimise water.
- Minimise drain blockages.

Quattro's unique design means a more efficient motor is used.


- Reduce energy use
- Reduce cost of repairs

Totals are calculated based on 30 cycles per day. "Washer" is a representative 1 pan bed pan washer machine. "Alternative" is a representative non-Haigh bed pan macerator machine.

## HAIGH: NATIONWIDE SUPPORT

- Training for care staff FREE
- Benefits vs Washer systems
- Options to suit your care environment
  - - \* Free user training and materials provided at same time as other on-site services purchased from Haigh
- Installation and Commissioning
- Minimising Total Lifetime Costs
- Genuine spare parts
- Consumables, Pulp and Chemicals
- Servicing, Maintenance and Call-outs
- Support for all macerator brands
- Asset surveys and management

01989 763131 info@haigh.co.uk

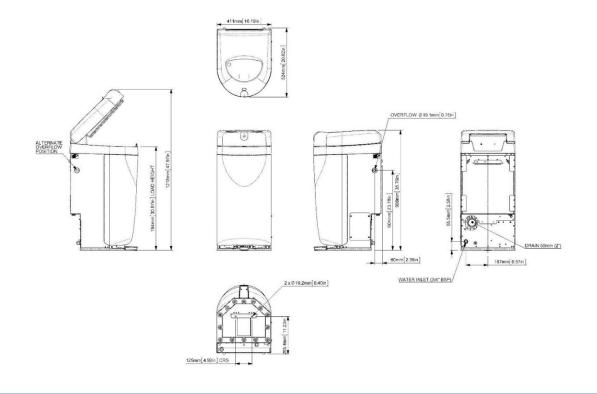

## **FAST – QUIET – EFFICIENT**

#### **Technical Performance**

Cycle performance: 82-180 seconds | 52-56 dBA | 54-59 dB max | Up to 4 standard medical pulp items

Hands free open | Option for hands free closing to improve infection control

Water: 15 - 22 litres per cycle used | 5.5 litres per minute cold water supply required

Electricity: < 0.01 kWh per cycle | 230V / 1Ph / 50Hz or 400V / 3Ph / 50Hz supply required

Inlet, Outlet & Cistern: 3/4" BSP Inlet | 2" BSP multi-fit P-trap Outlet | 11 litre cistern: EN 13077

Product Footprint: Width: 411 mm | Loading Height: 784mm | Depth: 524mm | Dry Weight: 58.5 kg

Haigh Engineering Ltd,
Alton Road,
Ross-on-Wye,
Herefordshire, HR9 5NG,
United Kingdom

Web: www.haigh.co.uk
Tel: 01989 763131
Email: info@haigh.co.uk















